




Technology Advancements in Hip Arthroscopy Visualization and Accessibility

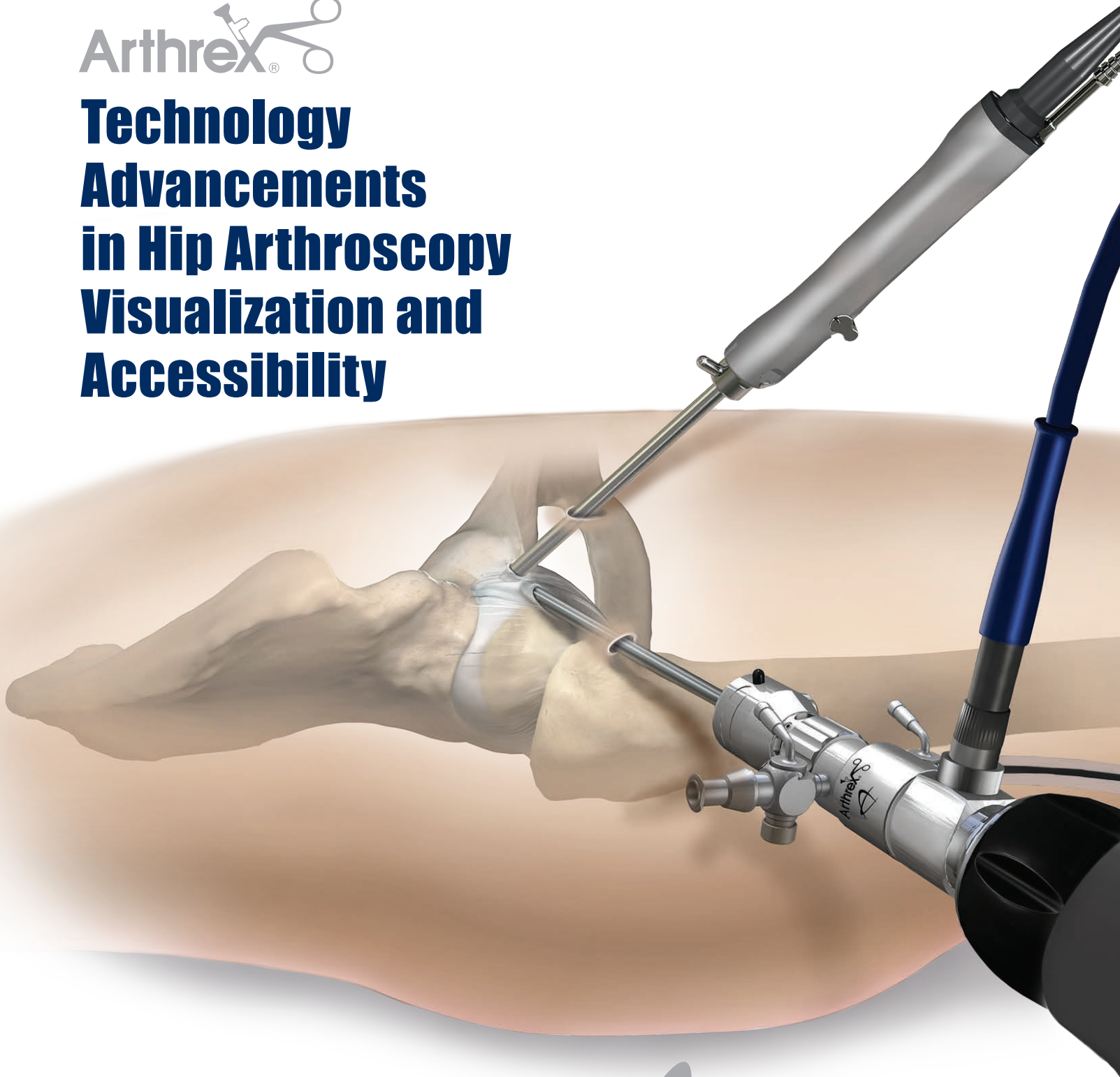


 **Synergy**^{UHD4™}

Optimal Hip Joint Access and Viewing Instruments

Advancements in Hip Joint Access


The Hip Length Arthroscope System was designed with a clear focus on accurate and safe hip joint access during arthroscopic procedures of the hip. It contains 30° and 70° arthroscopes in either a C-Mount or Eye Piece version with cannulated obturators that fit easily over Nitinol guide wires to minimize the insertion steps into the joint. The intuitive quick connect bridge allows for fast and simple portal switching between the atraumatic color-coded metal cannulas and the extended working length of the arthroscopes accommodates varying patient size, making this the most versatile hip arthroscope set available.

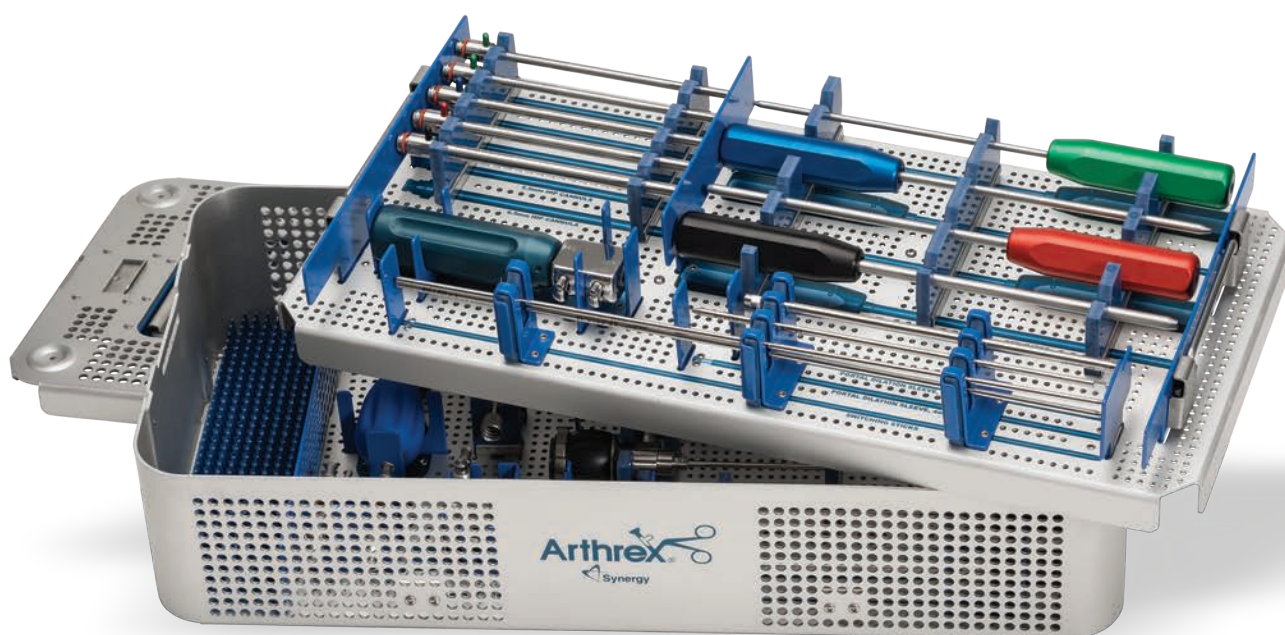

Synergy Console

30° and 70° Arthroscopes

The C-Mount and Eye Piece Hip Arthroscopes have been designed to meet the demands of varying patient anatomy. The high flow sheaths* are available in a single or double port version and have an accompanying cannulated obturator** to minimize the required steps for insertion into the joint.

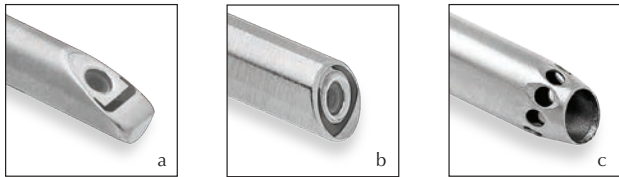


* Use with 3.5 mm Switching Sticks Only
 ** Use with either 1.1 mm or 1.5 mm Nitinol Wires




Atraumatic Metal Cannulas and Obturators

The nonfenestrated, atraumatic color-coded metal cannulas have a blunt tip to facilitate safe access into the joint space. The matching color-coded obturators have proximal threads for the metal cannulas to lock into for easier insertion through the soft tissues. The cannulas are available in half millimeter increments ranging from 4.5 mm - 6.5 mm diameters and have a working length of 145 mm.

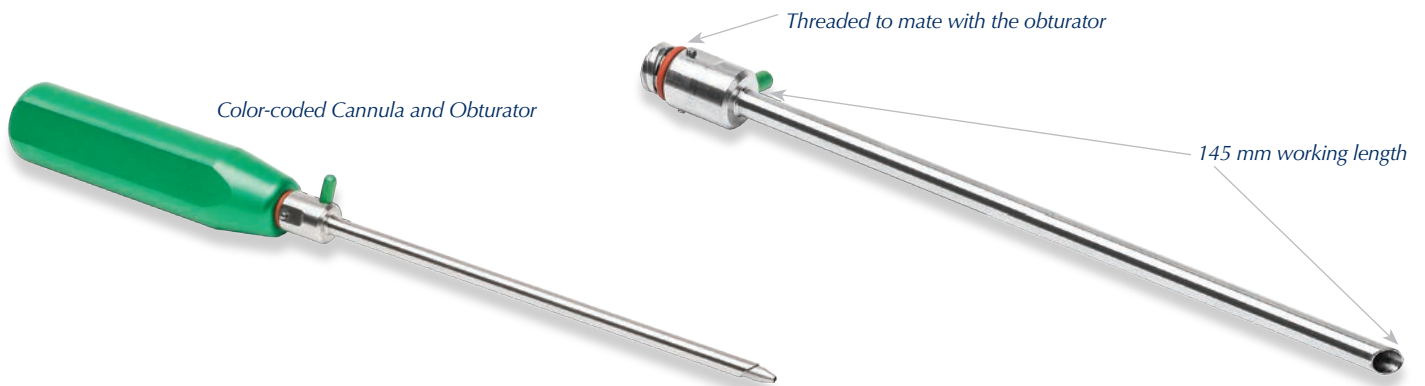

Quick Connect Hip Bridge System

The durable Quick Connect Hip Bridge is designed for fast and simple portal switching via a button connection to both the scope and atraumatic metal cannula. Available with either a single or double port, this bridge component provides an easy and reliable mechanism for fast and simple portal switching.

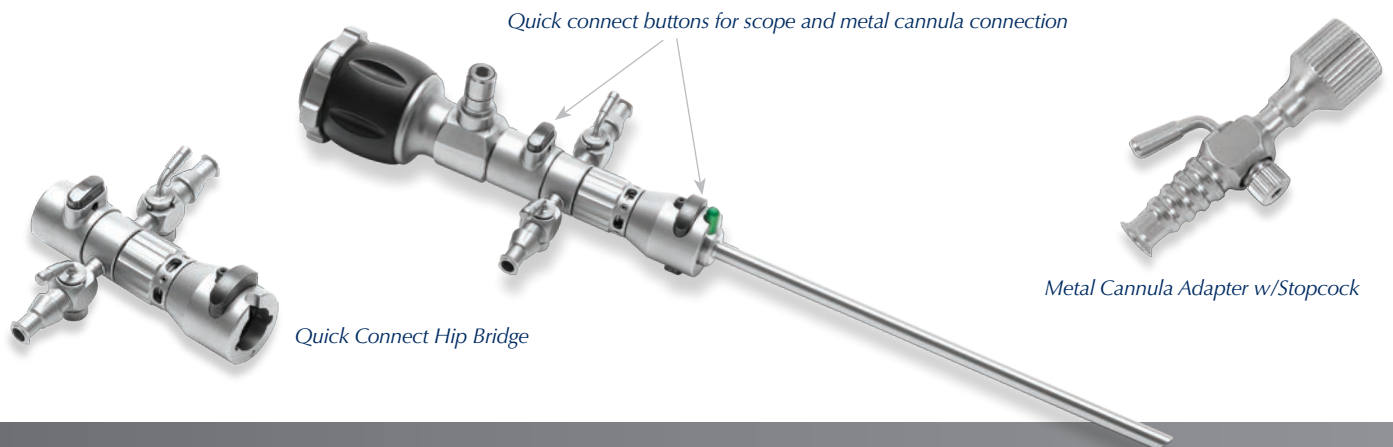


Joint Access Accessories

The sequential telescoping Portal Dilation Handle and Sleeves allow for fast and simple portal expansion to accommodate larger instruments and working cannulas. The tapered and blunt tip detachable Open Metal Cannulas can be removed from the Tear Drop Handle after insertion, minimizing instrument bulk outside the body.



Portal Dilation Handle and Sleeves

Open Metal Cannulas and Tear Drop Handle

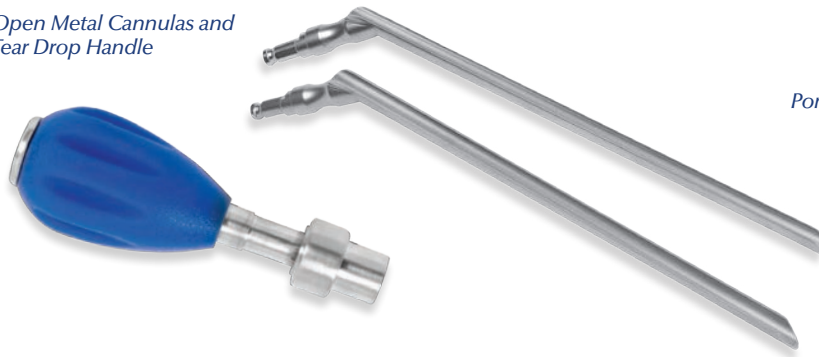

Hip Access Disposables

Arthrex offers three disposable kit options containing all the necessary components for precise hip joint access, as well as various blade options designed to provide a safe and controlled hip joint capsulotomy.

AR-6526-01S - Hip Access Disposable Kit



AR-6527S - Hip Arthroscopy Disposable Kit w/Banana Blade



AR-6526S - Hip Arthroscopy Master Disposable Kit w/Knife Blades

Retractable Cannulated Knife

Banana and Hook Blades



Ordering Information

Extended Length Arthroscope System AR-6576S-E (includes Eye Piece Scopes):

Hip Length Metal Cannula, 4.5 mm, nonfenestrated, qty. 2	AR-3420H-4.5
Hip Length Metal Cannula, 5.0 mm, nonfenestrated	AR-3420H-5.0
Hip Length Metal Cannula, 5.5 mm, nonfenestrated	AR-3420H-5.5
Hip Length Metal Cannula, 6.5 mm, nonfenestrated	AR-3420H-6.5
Hip Length Obturator, cannulated, 4.5 mm	AR-3421H-4.5
Hip Length Obturator, cannulated, 5.0 mm	AR-3421H-5.0
Hip Length Obturator, cannulated, 5.5 mm	AR-3421H-5.5
Hip Length Obturator, cannulated, 6.5 mm	AR-3421H-6.5
Switching Stick, cannulated, 4.0 mm x 15"	AR-6520-15
Switching Stick, cannulated, 3.5 mm x 15", qty. 2	AR-6520-16
Metal Open Cannula, tapered	AR-6507-03
Metal Open Cannula, blunt tip	AR-6507-04
Tear Drop Handle	AR-2001
Portal Dilation Sleeve, 4 mm	AR-6524-04
Portal Dilation Sleeve, 6 mm	AR-6524-06
Dilation Handle	AR-6524
Hi-Flow Sheath w/Double Stopcock, fenestrated, qty. 2	AR-3373-4002H
Conical Obturator for 3.5 mm Hip Scope Sheath, qty 2	AR-3375-4001H
Cannulated Obturator for 3.5 mm Hip Scope Sheath	AR-3375-4003H
Quick Connect Hip Bridge w/Double Stopcock	AR-3380-4004H
Inflow Cannula Adapter w/Stopcock	AR-3035L
30° Eye Piece Hip Arthroscope, 3.5 mm x 202 mm	AR-3350-4030H
70° Eye Piece Hip Arthroscope, 3.5 mm x 204 mm	AR-3350-4070H
Instrument Case	AR-6576C

Extended Length Arthroscope System AR-6576S-C (includes C-Mount Scopes):

Hip Length Metal Cannula, 4.5 mm, nonfenestrated, qty. 2	AR-3420H-4.5
Hip Length Metal Cannula, 5.0 mm, nonfenestrated	AR-3420H-5.0
Hip Length Metal Cannula, 5.5 mm, nonfenestrated	AR-3420H-5.5
Hip Length Metal Cannula, 6.5 mm, nonfenestrated	AR-3420H-6.5
Hip Length Obturator, cannulated, 4.5 mm	AR-3421H-4.5
Hip Length Obturator, cannulated, 5.0 mm	AR-3421H-5.0
Hip Length Obturator, cannulated, 5.5 mm	AR-3421H-5.5
Hip Length Obturator, cannulated, 6.5 mm	AR-3421H-6.5
Switching Stick, cannulated, 4.0 mm x 15"	AR-6520-15
Switching Stick, cannulated, 3.5 mm x 15", qty. 2	AR-6520-16
Metal Open Cannula, tapered	AR-6507-03
Metal Open Cannula, blunt tip	AR-6507-04
Tear Drop Handle	AR-2001
Portal Dilation Sleeve, 4 mm	AR-6524-04
Portal Dilation Sleeve, 6 mm	AR-6524-06
Dilation Handle	AR-6524
Hi-Flow Sheath w/Double Stopcock, fenestrated, qty. 2	AR-3373-4002H
Conical Obturator for 3.5 mm Hip Scope Sheath, qty 2	AR-3375-4001H
Cannulated Obturator for 3.5 mm Hip Scope Sheath	AR-3375-4003H
Quick Connect Hip Bridge w/Double Stopcock	AR-3380-4004H
Inflow Cannula Adapter w/Stopcock	AR-3035L
30° C-Mount Hip Arthroscope, 3.5 mm x 202 mm	AR-3355-4030H
70° C-Mount Hip Arthroscope, 3.5 mm x 204 mm	AR-3355-4070H
Instrument Case	AR-6576C

Ordering information continued on the back

Items Sold Separately:

Single Stopcock Hi-Flow Sheath	AR-3373-4001H
Quick Connect Hip Bridge w/Single Stopcock	AR-3380-4003H
Switching Stick, noncannulated, 3.5 mm x 15"	AR-3020-03
Switching Stick, noncannulated, long, 3.5 mm x 18"	AR-3020-04
Low Profile Conical Obturator for 3.5 mm Hip Scope Sheath	AR-3375-4007H

Hip Access Disposables:

Hip Access Kit w/18" 1.1 mm Nitinol Wires	AR-6526-01S
Hip Access Kit w/15" 1.5 mm Nitinol Wires	AR-6526S
Hip Access Kit w/15" 1.1 mm Nitinol Wires and Banana Blade	AR-6527S
Retractable Cannulated Knife	AR-6527-01
Banana Blade	AR-6527-05
Hook Blade	AR-6527-19



View U.S. patent information at www.arthrex.com/corporate/virtual-patent-marking

©2016, Arthrex Inc. All rights reserved. LB1-0027-EN_B