



PivotPost

Adjustable Movement Support

The PivotPost is a lateral post that moves in 3 planes to easily place it in the proper position in relation to the patient and pivot out of the way when not in use. The unique over-the-top sweeping design prevents the patient's leg from moving over the post when valgus stress or a figure 4 is applied.

Features and Benefits

PivotPost Movement

- Prevent patient's leg from moving over post when valgus stress is applied
- 3 planes of motion, including up/down **(a)**, in/out **(b)**, and rotation **(c)**
- Small size allows surgeons to operate easily; minimal weight allows for simple setup for OR staff
- Disposable pad that is easy to put on and take off
- Ergonomic, oversized buttons are conveniently located under the drape

Clark Rail Adapter **(d)**



- Custom-made of stainless steel for the PivotPost
- Finger tight mechanism holds firmly
- Pull knob to raise and lower is round; tightening knob is scalloped

Foam Pad **(e)**

- Five disposable (nonsterile) pads
- No specific top or bottom; feature a double-sided hook-and-loop fastener

Ordering Information

Product description	Item number
PivotPost	AR-1560
Knee post	AR-1560-02
Clark rail adapter	AR-1560-03
Foam pads (box of five)	AR-1560-01

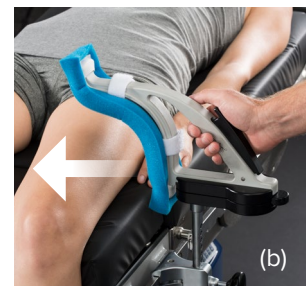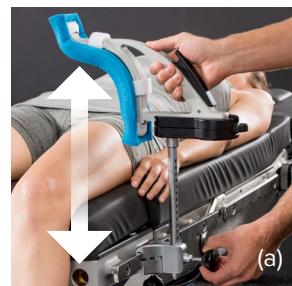

arthrex.com

© 2024-08 Arthrex, Inc. All rights reserved. LS1-0165-EN_C

